

Dansac NovaLife Features and Benefits

Closed Pouches 1 and 2 piece

Release tab

The barrier can be simply and securely removed thanks to the easy-to-see orange non-adhesive release tab.

Off centre starter hole

The off centre starter hole on the cutting guide means that the pouch sits lower on the abdomen, making it easier to hide.

Water repellent cover



The soft fabric of the pouch is water repellent and will stay dry even when bathing.

Oval shaped barrier

The unique oval shaped design suits the natural curves of the human body and ensures that the adhesive area is positioned where it is needed most.

Highly effective NovaLife filter

Membranes on both sides of the NovaLife filter protect it from the contents of the pouch as well as from water outside and make NovaLife filter patches unnecessary. The filter can expel a large amount of gas very quickly, which minimises the risk of the pouch ballooning.

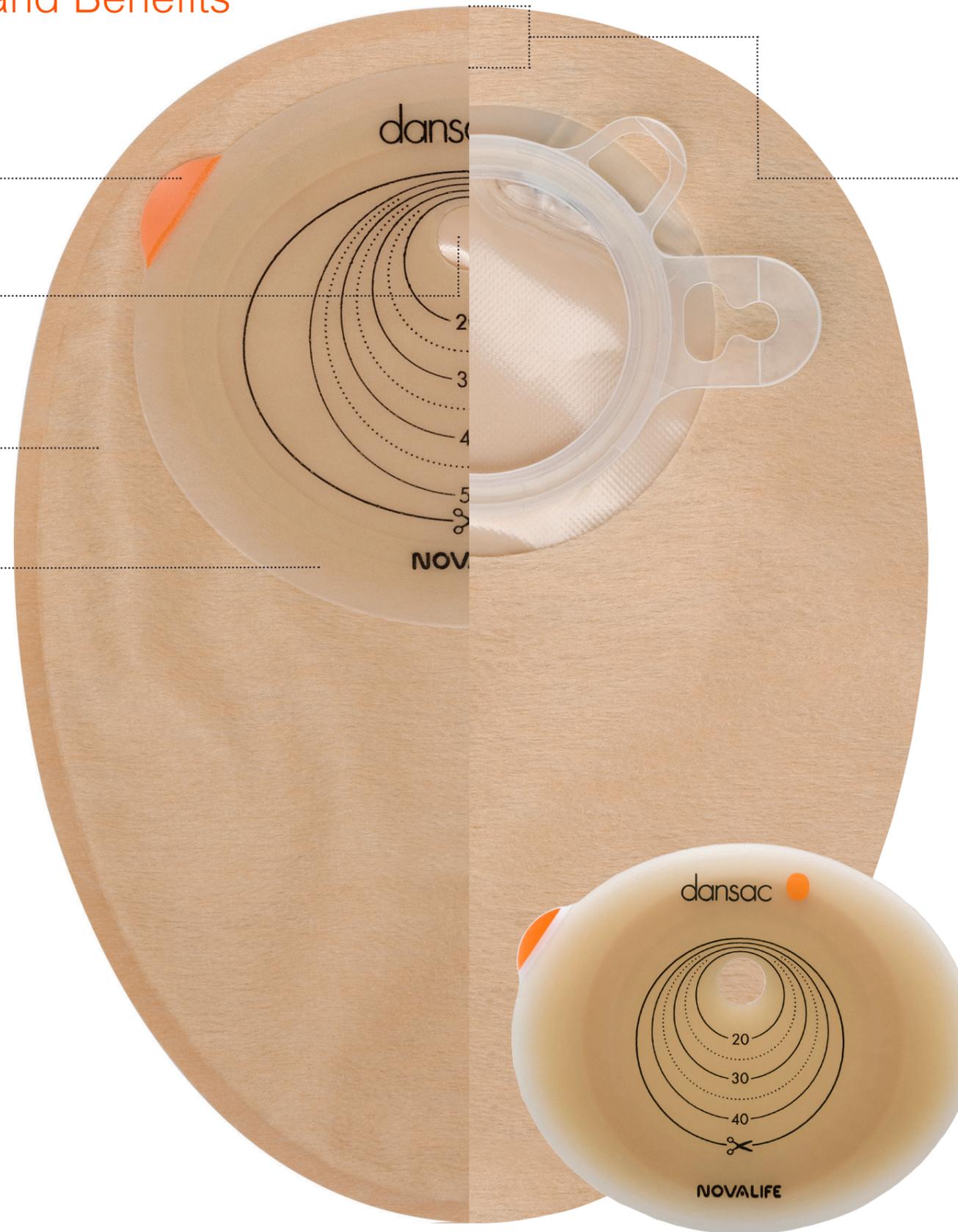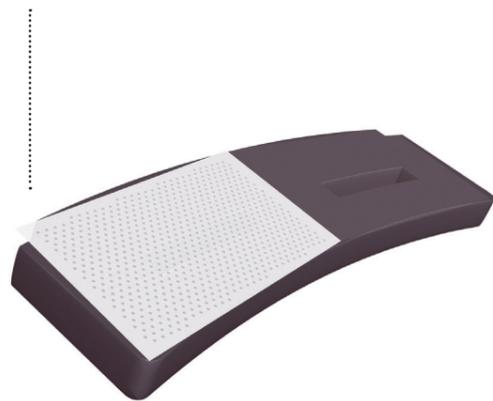

Low headspace

The pouch area above the barrier has been minimised, which reduces the height of the pouch thus makes it more discreet and prevents it from slouching.

Secure adhesion

The surface of the barrier is thickest around the stoma then thins towards the edges to provide secure and gentle adhesion to the skin.



Skin-friendly barrier

The highly skin-friendly composition of the GX Hydrocolloid barrier is based on natural ingredients and is proven over many years.

Secure and flexible coupling

The pouch and the wafer are securely connected by the proven coupling system. The ring moves freely to enhance the way the wafer follows the contours of the body.

